



# Cheminert UHPLC Valves and Selectors

## Notice of Design Change



As of early 2015, C7\* UHPLC injectors, switching valves, and selectors with 1/32" and 1/16" fittings have been upgraded to the C8\* design. The only modification is a longer body to accommodate a slightly different spring arrangement. The new arrangement provides more consistent spring force, resulting in better performance and lifetime.

The temperature and pressure specifications, port sizes, rotors, and stators remain identical.

Valves with 360 micron fittings remain unchanged, and retain the C7\* part number designation.

### Length dimensions (L)

Number of ports or positions	C7* (older)		C8* (current)	
	mm	in	mm	in
4	38.15 mm	1.5"	44.45 mm	1.75"
6	38.15 mm	1.5"	44.45 mm	1.75"
8	40.15 mm	1.58"	44.45 mm	1.75"
10	40.15 mm	1.58"	44.45 mm	1.75"

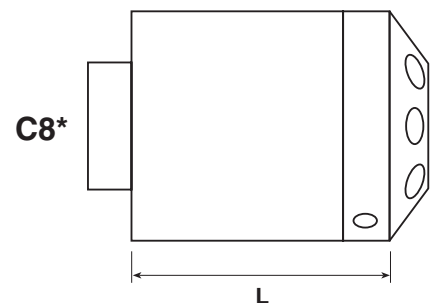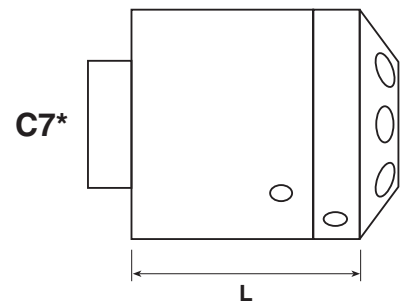

North America, South America, and Australia/Oceania contact:

**VICI** Valco Instruments Co. Inc.  
 P.O. Box 55603  
 Houston, TX 77255  
 Sales: (800) 367-8424  
 Tech: (713) 688-9345  
 Fax: (713) 688-8106 valco@vici.com

Europe, Asia, and Africa contact:

**VICI** VICI AG International  
 Parkstrasse 2  
 CH-6214 Schenkon  
 Switzerland  
 Phone: +41 41 925 6200  
 Fax: +41 41 925 6201 info@vici.ch  
 Cheminert® and VICI® are registered trademarks of Valco Instruments Co. Inc. and VICI AG