




Schematic Flow Diagram

Cheminert® Model C35ZT Trapping stream selectors

Cheminert Model C35ZT valves are used for multi-column, multi-sample, or multi-trap applications. A typical application, shown below, is the collection of fractions for analysis at a later time.

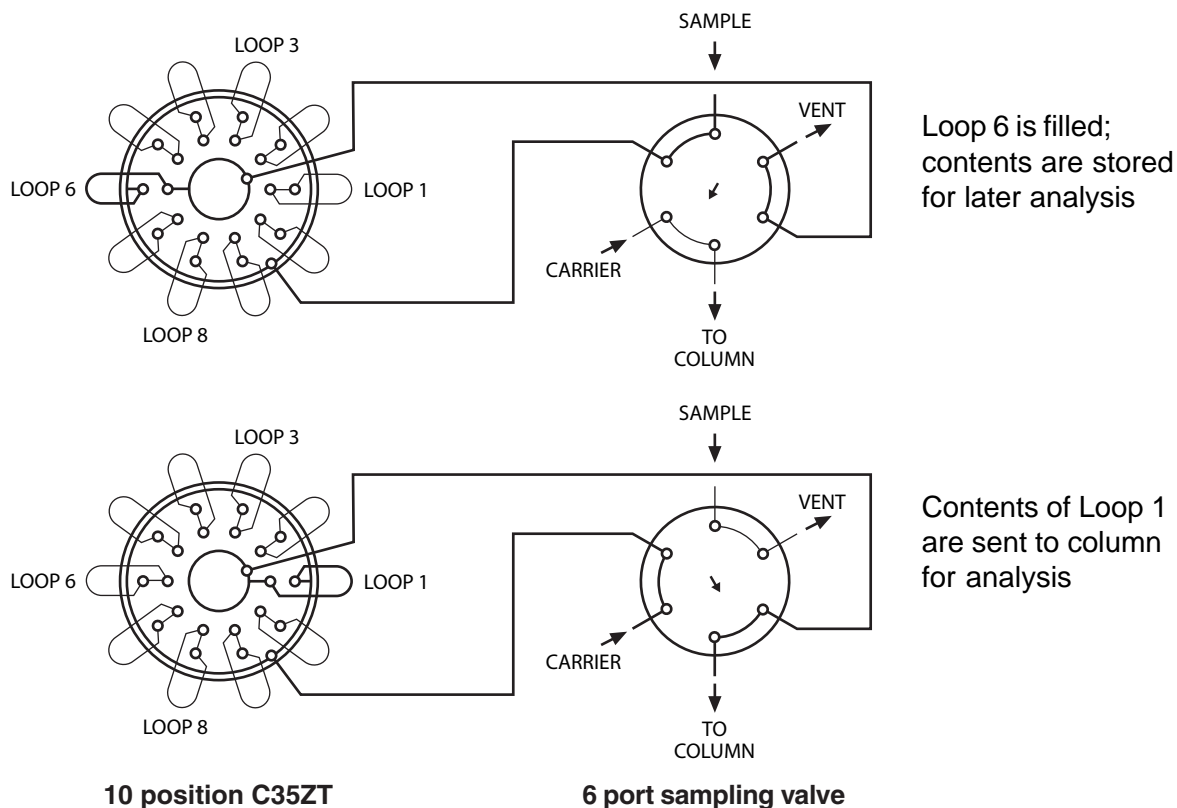

CAUTION:

The number of positions indexed by the actuator is *half* the quantity of numbered ports on the cap:

20 numbered ports + 1 outlet + 1 inlet = 10 positions

24 numbered ports + 1 outlet + 1 inlet = 12 positions

26 numbered ports + 1 outlet + 1 inlet = 13 positions

North America, South America, and Australia/Oceania contact:

VICI® Valco Instruments Co. Inc.
 P.O. Box 55603
 Houston, TX 77255
 Sales: (800) 367-8424
 Tech: (713) 688-9345
 Fax: (713) 688-8106 valco@vici.com

Europe, Asia, and Africa contact:

VICI® VICI AG International
 Parkstrasse 2
 CH-6214 Schenkon
 Switzerland
 Phone: +41 41 925 6200
 Fax: +41 41 925 6201 info@vici.ch

TN-823 9/09